Zapper K100

GEBRAUCHSANWEISUNG
OPERATING INSTRUCTIONS
ISTRUZIONI D'USO

DE	Gebrauchsanweisung	4
EN	Operating instructions	11
IT	Istruzioni d'uso	17

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Important Safety Instructions

Please read the operating instructions carefully before using the device!

- The Zapper K100 is designed for medium frequency treatments. Do not use the device for any other purpose.
- Only use the device with original accessories.
- Keep the device away from water and other liquids.
- Do not drop the Zapper K100, do not use it inappropriately or expose it to extreme temperatures or high humidity (only use at temperatures between 10 °C and 40 °C and a relative humidity below 90 %).
- Do not use the Zapper K100 if it is not working properly or if it has been damaged in any way.
- Always store the Zapper K100 in the original case after use to protect it from damage or dust.

Safety Guidelines

Patients with an implanted medical device (f. e. pacemakers) may only be treated with the Zapper K100, after consulting their physician first. The Zapper K100 may only be connected to one patient at a time.

Warnings!

The use of a high frequency surgery device at the same time can cause burns underneath the electrodes. The application of the device close to a short- or microwave device (~ 1 m) can cause variations in the output values of the Zapper K100.

Portable or mobile high frequency communication devices may have influence on electrical medical devices. Please read „Inauguration Of The Device“.

The electrodes shall not be smaller than 2 cm².

Description Of The Device

Upper surface:

1 Mode select switch


Front Panel:

2 Indicator lamp „Continuous Duty“

3-7 Indicator lamp „Program“

8 Exit hub

9 Indicator lamp

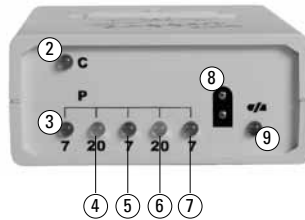Explanation Of Symbols

Attention: Read accompanying documents, especially user manual!

Type BF applied part

Date of manufacture

Catalogue number

Serial number

Environment: Do not dispose of the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this you will help to protect the environment.

CE0482 By labelling with CE mark, the manufacturer asserts that the product fulfills all relevant requirements of the respective EC Directives. A conformation assessment process has been successfully completed. The identification number in respect of the conformation assessment process is given in accordance with the CE labelling of the relevant notified body.

Technical Data

1-Channel medium frequency stimulator

Output current with 1 k Ω	8,7 V \pm 10 %
Frequency	32 kHz \pm 1 %
Pulse form	positive square
Pulse duty factor	1:1
Current consumption	\sim 7 mA
Power supply	9 Volt compound battery
Dimensions	\sim 115 mm x 60 mm x 30 mm
Weight (without accessories)	\sim 90 g
Classification	The Zapper K100 is classified in Class IIa according to directive 93/42/EEC, annex IX

Operating Of The Device

Before turning on the device, connect the battery (see „Battery Change“) properly, connect cables and electrodes and securely place electrodes on treatment area.

Turn on the device by sliding the mode select switch (1) to „P“ or „C“.

Please make sure that the Zapper K100 is not placed next to or

onto other electronic devices (f. e. mobile phones) in order to avoid reciprocal influencing. If this is impossible to avoid, pay attention to the functions while using the device and make sure the device is working properly.

Electrode Placement

Connect one electrode to each, left and right wrist. Connect the electrodes to the plugs on the end of each cable, then insert the double plug at the other end of the cables into the exit hub (8). If the electrodes are felt to be uncomfortable on the wrist, they can also be placed on the ball-and-socket footing. The polarity of the electrodes is not important.

In rare cases sensible patients may feel a prickling sensation underneath the electrodes. In this case a wet paper tissue can be placed between electrode and skin.

During the stimulation hands respectively feet shall not touch.

Set Mode

If the mode select switch (1) is set on „C“, a continuous stimulation without time limit is started. The indicator lamp „C“ (2) is blinking and an acoustic signal occurs at the beginning of the stimulation. If the mode select switch (1) is set on „P“ the following program will be started automatically:

7 Min. Stimulation – 20 Min. Pause – 7 Min. Stimulation – 20 Min. Pause – 7 Min. Stimulation

The indicator lamps (3), (4), (5), (6) and (7) will blink according to the program run. When starting or changing the program, an acoustic signal will occur.

Attention! The device does not turn off automatically after the end of the program. As long as the device will not be turned of, an acoustic signal will be given.

Electrode Control

The Zapper K100 is equipped with an automatic electrode control. If the electrodes are not connected properly or do not have enough contact to the patients skin, an acoustic warning signal will be given.

In this case, please check if the electrodes are connected and placed properly. If the warning signal still occurs, put a little bit of water on the skin or tissue underneath the electrode.

Battery Change

TO REPLACE THE BATTERY, TURN OFF THE DEVICE!

If the red battery light (9) flashes, the battery needs to be changed.

Slide off the cover of the battery compartment on the back of the device and disconnect the battery.

Check polarity of the battery when connecting the new 9 Volt compound battery (f. e. Type 6LR61) to ensure it is correctly attached.

If the device is not used for a prolonged period, remove the battery.

Where rechargeable batteries are used, please follow the instructions of the charging unit.

Dispose of old batteries properly.

Re-Settings, Modifications And Repairs

The manufacturer accepts responsibility for the safety and performance of the Zapper K100 provided that re-settings, modifications and repairs are carried out by persons authorised by him, and that the Zapper K100 is used in accordance with the operating instructions.

Warranty

We guarantee the Zapper K100, excluding cables, electrodes etc. for one year from the date of purchase.

Maintenance And Cleaning

No special detergents are necessary to clean the Zapper K100. The device should be cleaned with a soft, dust-free cloth.

Do not get the unit wet or expose to extreme humidity. However, if this occurs it must undergo a safety check before being used again.

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P GmbH
Pierenkemper GmbH